

Revolutionize your teeth with the **gemini®** Teeth Whitening System, in just **ONE hour!**
BEYOND® International Inc., global leaders in teeth whitening products





gemini®

Teeth Whitening Accelerator
Home Edition



- Next generation LED for maximum acceleration.
- No sensitivity.
- Removes years of stains caused by coffee, red wine and smoking. 6 to 8 shades Improvements!
- Embedded magnets line light to tray for a perfect fit every time.
- Simply Rechargeable via USB or outlet.
- Clinically proven to Kill germs.
- Refill kits available for continued maintenance!

Q&A

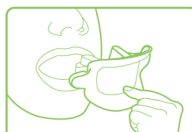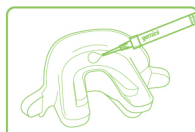

before after

What results should I expect from the Gemini® Teeth Whitening Accelerator?

Results will vary depending on tooth condition, but one should see on average a 6-8 shade improvement after one sitting of sixty minutes in total.

Can I expect tooth sensitivity?

No. Tooth sensitivity can occur from either high concentrations of hydrogen peroxide or dehydration from the gel or light. The Gemini® uses a mild 6% concentration, coupled with a UV-free LED light accelerator.

Do I have to do all 3-20 minute procedures in a row?

No. This is recommended for best results, but the procedure can be broken down to 20 minute sessions a day. We do recommend the full three sessions of twenty minutes in one sitting for the fastest and longest-lasting results.

What type of gel is being used?

We use a 6% hydrogen peroxide mix, which allows for a mild formula and yet provides strong results when combined with our specially designed Gemini® accelerator.

How long does the gel last?

The gel should last through the expiry date (around 18 months), but be sure to refrigerate the gel once it's opened.

Will the Gemini® work on porcelain veneers, bonded teeth and dental implants?

Veneers, bonding and dental implants all will collect surface stains over time. The Gemini® treatment will whiten these stains, but the whitening will not be as deep as with live, natural teeth. The whitening will not go past the initial color that was created for the material.

If I experience sensitivity, what are the next steps?

You should get in contact with your local dental professional to ensure your teeth and gums are healthy. It is unusual to have sensitivity with the Gemini® treatment, and this could mean you have another dental concern.

Can I use other gels with the Gemini® system?

The Gemini® gels were specifically formulated to safely react with the Gemini® accelerator. Extensive testing was used to provide premium whitening results without sensitivity. We cannot offer an opinion on other gels regarding safety or effectiveness.

I followed the steps, and would like for my sister to try also. Can I buy trays separately?

Yes, we offer both the tray and refill gels separately. You can contact us via our website at www.beyonddent.com to reorder.

Is the Gemini® FDA approved?

The Gemini® is listed as a Class I medical device with the FDA.

beyond® International Inc. — *Global leaders in aesthetic dental products.* www.beyonddent.com

BEYOND USA Headquarters
711 Julie Rivers Dr.
Sugar Land, TX 77478 USA
TEL +1-281-277-4352
E-mail: info@beyonddent.com

BEYOND European Headquarters
Piotra Skargi 12
65-416 Zielona Gora, Poland
TEL: +48 512 401057
Email: europe@beyonddent.com

BEYOND Asia Pacific Headquarters
Chaoyang District, Beijingggg
China 100029
E-mail: asia@beyonddent.com

Questions? Call +1-855-956-6668 M-F 9:00 a.m. – 5:30 p.m. CST

