





Teeth Whitening Accelerator

Home Edition

BY-GM224

-Instruction Manual-

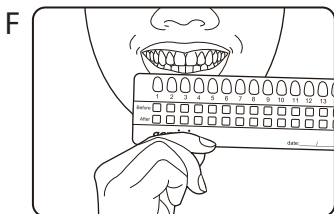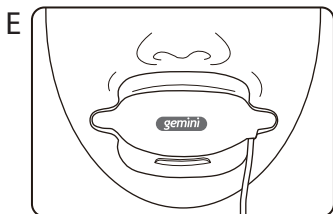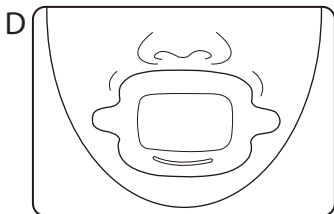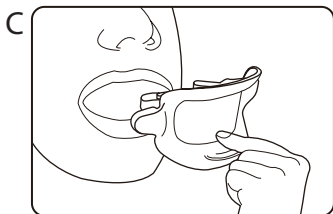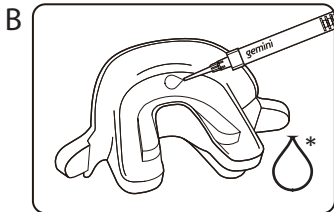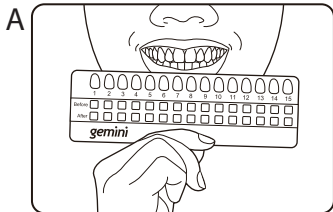English	3
Čeština	7
Deutsch	11
ελληνικά	16
Español	21
Français	26
Hrvatski	31
Italiano	36
Magyar	41
Nederlands	46
Norsk	51
Polski	55
Português	59
Romana	64
Русский	69
Slovenčina (Slovenský jazyk)	74
Svenska	79
Türkçe	84
简体中文	89
繁體中文	93
日本語	97
ภาษาไทย	102
한국의	106
ةيبرعلا	111

Check Your Package To Ensure That The Following Items Are All Included:

- 1 x Gemini™ LED light unit
- 5 x Gemini™ whitening gel syringes
- 1 x Gemini™ tray
- 1 x Gemini™ power unit with rechargeable lithium battery
- 1 x Gemini™ charge cord
- 1 x Gemini™ USB/electrical adaptor
- 1 x Bleaching Shade Guide

Before Use

Gemini™ light unit will require a full charge before initial use. For ease of use, the Gemini™ comes with a USB charging cord with adaptor. Via electrical outlet it will take approximately 3 hours to obtain a full charge, and can perform approximately twenty 20-minute cycles from a full charge. Charging via computer will be slightly longer, depending on system.

Light Indication	Status
Blue	Fully charged
Blue, pulsing slowly	In use
Orange	Charging
Orange, pulsing quickly	Needs recharging

Gemini™ Instructions For Use

1. Floss and brush teeth gently.
2. Determine original tooth shade using the Bleaching Shade Guide. (Fig. A) Mark shades for post-treatment comparison.
3. Using one of the included gel syringes, apply the whitening gel evenly along the margin of the tray between the inner ridge and the outer plastic lining, according to illustration. User should expend 50% of syringe gel to lower part of tray, then flip tray and expend remaining gel to upper part of tray. (Fig. B)
* Actual size dose
4. Insert tray inside mouth. Horizontal tab should be facing from the bottom of the tray. (Fig. C)

Bite gently. Teeth should fit over teeth outline on tray. Adjust the lips up and away from the tray and rest the lips in the lip guard portion of the tray for a comfortable fit. (Fig. D)

5. Place the Gemini™ light emitter in front of tray. Embedded magnets will automatically direct the light to the optimal position in front of the tray. (Fig. E)

Press and hold “” button for five seconds. Gemini™ light will activate for a full twenty minute cycle. At end of session, light will automatically deactivate.

If session needs to end prematurely, press and hold “” button for two seconds. Light will deactivate.

After twenty minute cycle, remove tray and rinse. Tray can be rinsed and patted dry with tissue or soft cloth.

Once dry repeat steps 3-5 two more times.

6. Evaluate color change against the marked shades from step 2. (Fig. F)

For Best Results

Minimize consumption of colored food and drinks while using the Gemini™. Try to limit use of tobacco products. Maintain good oral hygiene during the treatment.

In rare cases, whitening gel may come into contact with the gums during the treatment. In such a case the pink tissue may turn white post-treatment. This will fade in a few hours following the treatment and is not cause for concern.

Keep out of reach of children. Not suitable for pregnant or nursing women or children under 18 years of age.

If you have severe stains caused by tetracycline, fluoride, or another chemical, discuss your best options for whitening these stains with your dentist. The Gemini™ is effective on whitening these types of stains; however additional treatment cycles and adaptations may be necessary to achieve best results.

To help maintain your great new smile you may try the following products.

More information on our oral care line can be seen here: www.beyonddent.com

Cleaning The Gemini™

1. Wash the tray with soapy water, then rinse and dry with tissue or soft cloth.
2. Use a damp cloth to wipe off the Gemini™ LED cold-light unit. DO NOT rinse or submerge.

Storage And Maintenance

Store all items in a cool, dry place. Avoid exposure to direct sunlight. For best results, keep Gemini™ whitening gel refrigerated. After use, empty syringes may be recycled.

Note: Within the US, you may provide BEYOND® Dental & Health with your empty syringes for a \$5 coupon on your next order of Gemini™ refill syringe kit. Please contact BEYOND® Dental & Health for further details.

Contraindications

Gemini™ Whitening Gel is contraindicated for patients with known allergy or chemical sensitivity to peroxides or glycols.

Frequently Asked Questions

What results should I expect from the Gemini™ Teeth Whitening Accelerator?

Results will vary depending on tooth condition, but one should see on average a 6-8 shade improvement after one sitting of sixty minutes in total.

Can I expect tooth sensitivity?

No. Tooth sensitivity can occur from either high concentrations of hydrogen peroxide or dehydration from the gel or light. The Gemini™ uses a mild 6% concentration, coupled with a UV-free LED light accelerator.

Who can use the Gemini™ Teeth Whitening Accelerator?

The Gemini™ is user-friendly enough for anyone to use, but is generally suggested for those 18 or over. We also do not suggest either pregnant or nursing women to use.

Do I have to do all 3-20 minute procedures in a row?

No. This is recommended for best results, but the procedure can be broken down to a single 20 minute session per day. We do recommend the full three sessions of twenty minutes in one sitting for the fastest and longest-lasting results.

Do results vary per person?

Yes, the results will depend on your teeth condition and personal habits. Generally speaking, the Gemini™ will bring your teeth shade back to where it was originally before you stained it from tobacco use, coffee, red wine or staining food.

How long should I wait between procedures?

After a full teeth whitening, we generally suggest you wait 3-6 months before another session. It will generally take this long to re-stain your teeth from normal life habits.

What type of gel is being used?

We use a 6% hydrogen peroxide mix, which allows for a mild formula and yet provides strong results when combined with our specially designed Gemini™ accelerator.

How long does the gel last?

The gel should last through the expiry date (around 18 months), but be sure to refrigerate the gel once it's opened.

Will the Gemini™ work on porcelain veneers, bonded teeth and dental implants?

Veneers, bonding and dental implants all will collect surface stains over time. The Gemini™ treatment will whiten these stains, but the whitening will not be as deep as with live, natural teeth. The whitening will not go past the initial color that was created for the material.

If I experience sensitivity, what are the next steps?

Stop Treatment. You should get in contact with your local dental professional to ensure your teeth and gums are healthy. It is unusual to have sensitivity with the Gemini™ treatment, and this could mean you have a separate dental concern.

Can I use other gels with the Gemini™ system?

The Gemini™ gels were specifically formulated to safely react with the Gemini™ accelerator. Extensive testing was used to provide premium whitening results without sensitivity. We cannot offer an opinion on other gels regarding safety or effectiveness.

I followed the steps, and would like for my sister to try also. Can I buy trays separately?

Yes, we offer both the tray and refill gels separately. You can contact us via our website at www.beyonddent.com to reorder.

Is the Gemini™ FDA approved?

The Gemini™ is listed as a Class I medical device with the FDA.

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